SMOOTH SKIN
bare+

IPL USER GUIDE
GUÍA DEL USUARIO

10 minutes and *voilà!*

CONTENTS

ENGLISH

Page 5

ESPAÑOLA

Página 30

SMOOTHSKIN
bare+

CONTENTS

1. INTRODUCTION	5	8. HOW TO USE THE SMOOTHSKIN BARE+	21
2. HOW IT WORKS	6	9. CLEANING, MAINTENANCE AND STORAGE	26
3. IMPORTANT SAFETY INFORMATION	7	10. TRAVELING WITH THE SMOOTHSKIN BARE+	26
4. POTENTIAL RISKS, SIDE EFFECTS AND SKIN REACTIONS	17	11. TROUBLESHOOTING	27
5. BENEFITS OF USING THE SMOOTHSKIN BARE+	19	12. WARRANTY	27
6. WHAT IS THE INTENDED USE OF THE SMOOTHSKIN BARE+	21	13. LABELS AND SYMBOLS	28
7. YOUR SUITABILITY	21	14. TECHNICAL SPECIFICATIONS	29

Welcome to SMOOTHSKIN BARE+

Congratulations!


You've opened up a world of unparalleled beauty, smoothness and convenience. Come and join in the conversation at the beauty room 'support programme' www.smoothskin.com where you can swap tips and get expert advice throughout your treatment.

Package Contents

Inside the packaging of the SmoothSkin Bare+ you will find the following:

- Handset and Power Supply
- Mains Power Cable
- User Manual


The SmoothSkin Bare+ System

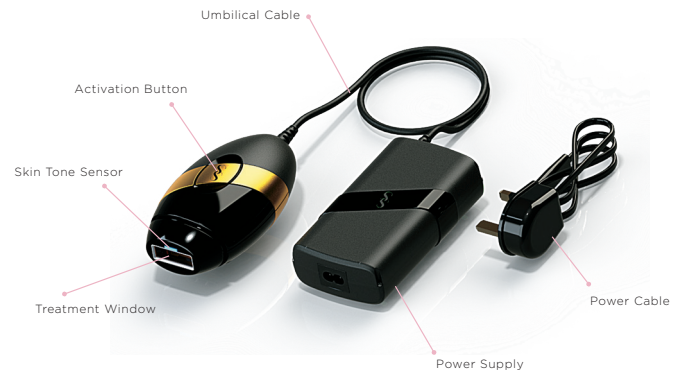2. HOW IT WORKS

HOW DOES IPL WORK?


SmoothSkin Bare+ is designed to help break the hair cycle. Light energy is transferred through the skin's surface and is absorbed by melanin present in the hair shaft. That light energy is converted to heat energy (below the surface of the skin), which stops the hairs from growing back. Now you can look forward to seeing the treated hairs naturally fall out over the course of a few days to 1-2 weeks.

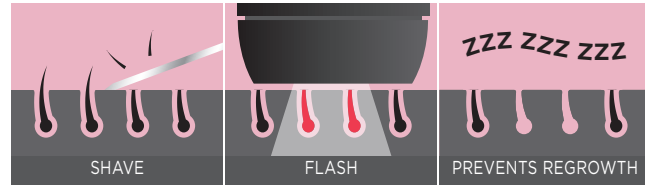3. IMPORTANT SAFETY INFORMATION

CONTRAINDICATIONS

DO NOT USE SmoothSkin Bare+ if you have very dark skin. The Skin Tone Sensor located near the output window of the device will determine the tone of your skin. If the Skin Tone Sensor determines that your skin is too dark for safe use, operation of the SmoothSkin Bare+ will be prevented.

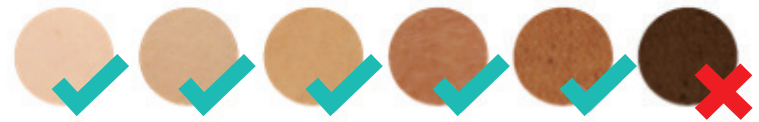DO NOT USE SmoothSkin Bare+ if your skin has been artificially or naturally tanned or is burned from sun overexposure. Your skin may be extra sensitive following sun exposure and particularly susceptible to burning, blistering, discolouring or scarring if you use the SmoothSkin Bare+. Wait 1 week until sunburn or tan disappears before using the SmoothSkin Bare+.

DO NOT expose treated areas to the sun. Wait at least 7 days after treatment before exposing treated skin to direct sunlight. Your skin may be extra sensitive following SmoothSkin Bare+ treatment and particularly susceptible to sunburn. Apply sunscreen (SPF 15 or greater) to the treated area or cover the treated area with suitable clothing.

DO NOT USE SmoothSkin Bare+ directly on the nipples, genitals or around the anus. These areas may have a darker skin colour and/or greater hair density. Using the device in these areas may cause discomfort/pain or injure (burn, discolour or scar) your skin.

DO NOT USE SmoothSkin Bare+ if you have a history of skin cancer or pre-cancerous lesions (e.g. nevi or a large number of moles)

DO NOT USE SmoothSkin Bare+ if you are pregnant, lactating, as this device has not been tested with these individuals.

DO NOT USE SmoothSkin Bare+ on dark brown or black spots such as birthmarks, moles or warts in the area you wish to treat. Using the SMOOTHKIN BARE may injure your skin or make existing conditions worse. You may experience side effects such as burns, blisters and skin colour changes or scarring.

DO NOT USE SmoothSkin Bare+ if you have any skin condition in the treatment area, including psoriasis, vitiligo, eczema, acne, herpes simplex or active infections or wounds.

DO NOT USE SmoothSkin Bare+ if you have a tattoo or permanent makeup in the area you wish to treat.

DO NOT USE SmoothSkin Bare+ if you are under the age of 18 as this device has not been tested on these individuals.

WARNINGS

EYE SAFETY

The SmoothSkin Bare+ emits flashes of Intense Pulsed Light. Direct exposure is potentially harmful to your eyes. Take care to follow the safety precautions below.

The Handset can only be activated (flashed) if the Skin Tone Sensor (next to the Treatment Window) detects a valid Skin Tone reading and is in full contact with the treatment area. However:

DO NOT look directly into the Treatment Window of the Handset when the system has electrical power and is switched on.

DO NOT attempt to activate (flash) the system towards the eye.

DO NOT treat areas around the eye (eyebrows or eyelashes).

ONLY attempt to activate (flash) the system when the front of the Handset is in good contact with the area you wish to treat.

WE RECOMMEND you look away from the Handset when firing onto your skin.

ELECTRICAL AND SAFETY

As with any electrical device, certain precautions must be taken in order to ensure your safety.

DO NOT try to open the SmoothSkin Bare+, as this may expose you to dangerous electrical components.

DO NOT USE if SmoothSkin Bare+ is damaged; for example, cracked Handset, cable damage (internal wires visible), cracked or broken glass on the Handset, etc.. This may result in injury.

DO NOT USE near bathtubs, showers, basins or other vessels containing water. If SmoothSkin Bare+ becomes wet, it should not be used.

DO NOT USE if the SmoothSkin Bare+ becomes uncomfortably hot to touch. This may indicate that the unit is damaged.

DO NOT USE a Power Supply or other component, other than those provided for use with SmoothSkin Bare+. This may damage the unit and cause it to stop working.

DO NOT attempt to activate (flash) the SmoothSkin Bare+ onto any surface other than skin.

WARNING no modification of this equipment is allowed.

DO NOT position SmoothSkin Bare+ so that it is difficult to disconnect from the Mains Supply.

GENERAL SAFETY

DO NOT USE if you have had a skin peel treatment on the area you wish to treat. Please wait 30 days after a skin peel before treating with the SmoothSkin Bare+.

Treating tanned skin with Intense Pulsed Light (IPL) may result in permanent Hyperpigmentation (Darkening) or Hypopigmentation (lightening) of the treated skin.

DO NOT USE If you have a known sensitivity to sunlight (photosensitivity) or are taking medication that makes the skin more sensitive and causes photosensitivity, e.g. Retin A, Accutane and/or other topical retinoids. Always check the instruction leaflet that comes with your medicinal

product, to see if photosensitivity is listed as a side effect. Use of the SmoothSkin Bare+ on photosensitive skin could result in skin damage such as swelling or blisters.

The SmoothSkin Bare+ should be kept out of the reach of children under 18 years of age. Children should not use, play with, maintain or clean the device. There are a number of potential risks – exposure to the light output causing burns or eye damage, electrocution, strangulation from the cables.

The SmoothSkin Bare+ should be kept out of reach of pets and pests which could damage the device. Potential risks include reduced function and electrocution.

The SmoothSkin Bare+ is not intended for use by persons with reduced physical, sensory or mental capabilities. Unless they have been given supervision or instruction concerning safe and proper use of the device

by a person responsible for their safety.

The SmoothSkin Bare+ Glass Filter is a critical component that will get hot during use.

DO NOT touch the Filter during or directly after use. Regularly inspect the filter for damage. **DO NOT** use SmoothSkin Bare+ if the filter is cracked or missing.

PRECAUTIONS

Hair removal by lasers or intense pulse light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage treated on the face and neck.

DO NOT USE the SmoothSkin Bare+ on any area where you may want hair to grow back.

DO NOT USE the SmoothSkin Bare+ device on man's face, jaw, or neck. The hair in a man's beard area is too dense, and using the SmoothSkin Bare+ device in these areas may result in skin injury.

SmoothSkin Bare+ is not recommended for use on red, grey or blonde hair. The SmoothSkin Bare+ is not effective on these hair colours.

4. POTENTIAL RISKS, SIDE EFFECTS AND SKIN REACTIONS

Use of the SmoothSkin Bare+ could result in side effects.

The table below lists the known skin reactions that may occur following treatment with the SmoothSkin Bare+.

SIDE EFFECT	HOW TO ASSESS & REACT
Mild pain / discomfort in the area being treated.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed, and the pain should diminish with continued use.
Warm feeling or tingling sensations during treatment, which typically disappear after a few seconds to a minute and decrease with continued use.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed, and the pain should diminish with continued use.
Itchiness in the treatment area.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed.
Skin redness during or after treatment which disappears within minutes.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed once skin redness has disappeared.
Skin redness which does not disappear within minutes and lasts 24-48 hours.	Stop using the device immediately and consult your physician before using it again.
Pain or discomfort that is intense during treatment or persists after a treatment.	Stop using the device immediately and consult your physician before using it again.

IN VERY RARE CASES	HOW TO ASSESS AND REACT
Swelling and redness that does not disappear within two to three days.	Stop using the device immediately and consult your physician before using it again.
Temporary changes in skin colour (lightening or darkening).	If your skin colour changes, stop using the device immediately and consult your physician.
Blistering or burning of the skin.	Stop using the device immediately and consult your physician before using it again. Apply a cold pack to the area. Treat with antiseptic or burn cream.

In a clinical trial of 50 subjects, each subject received 12 weekly treatments. The recorded effects were:

- **Itchiness** – Four subjects experienced mild itchiness after treatment.
- **Slight Redness** – This was reported by 1 subject.
- **Mild Stinging Sensation** – Five subjects reported mild discomfort following treatment.
- **In-growing Hair** – One subject reported experiencing an in-growing hair following treatment.
- **Follicle Redness** – One subject reported a red spot at the treatment site following treatment.

In all cases, the treated skin returned to normal within 7 days.

5. BENEFITS OF USING THE SMOOTHSKIN BARE+

The SmoothSkin Bare+ is indicated for the permanent reduction of unwanted hair.

A clinical trial was performed to assess both the safety and efficacy of the SmoothSkin Bare+. The key features of the clinical trial were:

- All treatments were performed by non-medically trained operators at a single location using devices identical to the SmoothSkin Bare+ provided in this package.
- All subjects were required to complete questionnaires related to their general health and those who could not be treated, e.g. those who had one or more of the contraindications as listed in Section 3 of this User Manual were excluded. In addition, all subjects were required to provide Informed Consent as per international clinical trial requirements.
- A total of 53 female subjects, aged between 19 and 45 years old, commenced the trial.
- Three subjects left the trial during the treatment period for personal reasons. No subjects left the trial due to problems associated with the use of the SmoothSkin Bare+.
- Each subject received 12 weekly treatments to their chosen body location. A typical treatment consisted of shaving the site then applying the SmoothSkin Bare+ following the same process as outlined in Section 7, Step 3 of this user manual.
- Hair counts were used to measure the change in the amount of hair 6 months after the last treatment and 12 months after the last treatment. The hairs in the treatment sites were counted via high resolution photography and the difference in hair counts was calculated as a percentage change.
- When used as directed, the clinical trial participants showed on average 44% less hair 6 months after their last treatment and 36% less hair 12 months after their last treatment when compared to hair counts taken before treatment (See Table below). Actual results did vary from person to person.

SmoothSkin Bare+	
Number of Subjects at 6 Months Post Treatment	50
Hair Reduction at 6 Months Post Treatment	43.9%
Number of Subjects at 12 Months Post Treatment	33
Hair Reduction at 12 Months Post Treatment	36.0%
% of Subjects met success (>30% hair reduction) on all body areas at 12 months post-treatment. Subject Success is defined as greater than 30% hair reduction at all treatment sites at 12 months.	66.7%

- The incidence of side effects during the clinical trial was minimal (as outlined in Section 4 above) and the majority of subjects (48 out of 50) would recommend the SmoothSkin Bare+ to their friends.

6. WHAT IS THE INTENDED USE OF THE SMOOTHSKIN BARE+

INTENDED USE

The SmoothSkin Bare+ is indicated for the removal of unwanted hair. The SmoothSkin Bare+ is also indicated for the permanent reduction in hair regrowth, defined as the long-term, stable reduction in the number of hairs re-growing when measured at 6, 9 and 12 months after the completion of a treatment regime.

7. YOUR SUITABILITY


BODY AREAS

SmoothSkin Bare+ is suitable for use on body hair and female facial hair below the cheekline.

SKIN COLOUR

The SmoothSkin Bare+ is suitable for use on light, medium and dark skin tones up to and including skin tone 5. The SmoothSkin Bare+ incorporates a Skin Tone Sensor that measures the colour of your skin and will prevent use if your skin tone is too dark.

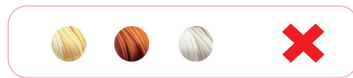HAIR COLOUR

The SmoothSkin Bare+ is suitable for use on naturally black or brown hair. It may not be as effective on white, grey, blonde or red hair.

8. HOW TO USE THE SMOOTHSKIN BARE+

WARNING: Follow each step in the treatment process ensuring you complete all aspects of each step before progressing to the next.

Step 1: PREPARE the area to be treated


Remove all visible hair by shaving. Make sure there is no hair remaining above the surface of the skin to ensure optimal treatment and prevent debris from covering the front of the SmoothSkin Bare+.

Cleanse the treatment area and pat dry. Remember, don't use any gels or creams with this device!

Step 2: Doing a TEST PATCH

Before your first treatment on each new body area, we RECOMMEND you test your skin in that area for a reaction to SmoothSkin Bare+.

The patch test area should be approximately 3cm x 2cm in size (equivalent to 2 flashes applied to the skin side by side). Follow Step 3 to treat the area.

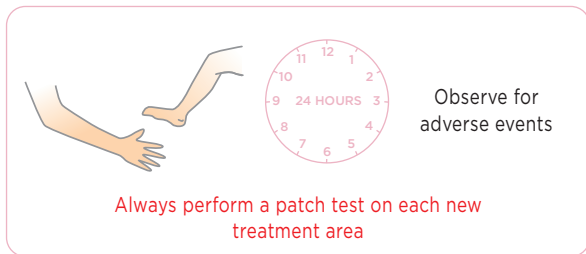Wait 24 hours following the patch test, to ensure your skin is suitable for treatment and there is no adverse reaction to the light energy. If there is no reaction after 24 hours, you may treat the area around the patch test. The 'patch tested area' should NOT be re-treated for at least 1 week.

Step 3: Treating with the SmoothSkin Bare+

A) Plug the Power Supply into an Electrical Outlet using the provided Power Cable. The Skin Tone Sensor on the front of the Handset will illuminate blue and the fan will start up.


STANDBY MODE Timeout

When the unit is left "inactive" i.e. no buttons are pushed for 10 minutes then the unit will revert to STANDBY mode, and the STS lights will switch off. To return to READY mode for treatment, press the Activation Button.

B) Press the Handset firmly against the area to be treated ensuring the Skin Tone Sensor is in contact with the skin. If the Skin Tone Sensor does not detect a valid skin tone, then the Handset cannot be activated. When the Skin Tone Sensor measures a valid skin tone, the lights around the Activation Button will illuminate WHITE to indicate the device is ready to flash.

- C) Press the ACTIVATION BUTTON. The device will flash, and a warm sensation should be felt.
- D) Move the Handset to the next area to be treated and repeat. The device needs to recharge between flashes. This will only take about 0.6 seconds.
- E) Ensure that the ENTIRE AREA receives treatment, and treat any stubborn areas (such as underarms) with multiple passes (running your device over the same area a maximum of 3 times).

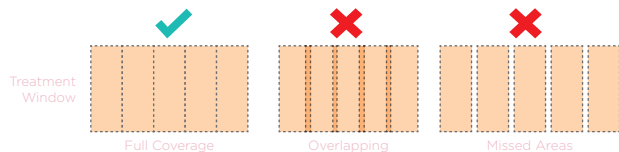WARNING: Make sure not to cover the vents on the device whilst treating, as this can lead to the device overheating.

SmoothSkin Bare+ is an extremely fast and powerful device, allowing you to complete a whole body treatment in under 10 minutes.

TREATMENT METHOD - STAMP OR GLIDE?

Stamp Method

This method is recommended for small areas such as the underarm and bikini line. Place the Handset on the skin, press & release the Activation Button. Remove from the skin before carefully re-applying to the next treatment area.

Glide method

This method is recommended for larger areas such as legs. Place the Handset on the skin, hold down the Activation Button (press down continuously), then glide the Handset along the skin between flashes. Try to move the Handset at a steady speed that gives full coverage without overlap or missed areas.

Step 4: Your SmoothSkin Bare+ Treatment Schedule

The SmoothSkin Bare+ should be used weekly for at least 12 weeks or until the desired results are achieved.

Step 5: Switching Off the SmoothSkin Bare+

When the treatment is finished, SmoothSkin Bare+ should be unplugged from the Electrical Outlet.

Refer to **section 9** for details of how best to look after your SmoothSkin Bare+.

9. CLEANING, MAINTENANCE AND STORAGE

After treatment, always switch off SmoothSkin Bare+ by unplugging from the Electrical Outlet.

After use, the Handset, Skin Tone Sensor and Treatment Window should be inspected for damage and wiped down with a dry, lint free cloth.


WARNING: The glass filter can get hot during use. **DO NOT** touch or clean the filter for at least 5 minutes after use to allow it to cool down.

NEVER use water or other cleaning fluids, as these can damage the device leading to a potential safety hazard.

If necessary, the Reflector and Skin Tone Sensor can be carefully cleaned by using a slightly damp cotton bud.

Store SmoothSkin Bare+ in a cool, dry place. Make sure that the Treatment Window and Skin Tone Sensor on the Handset are protected from damage.

Regularly check the device (including cords) for visible signs of damage. In case of damage or cracks, stop using the device - and visit the website for support details.

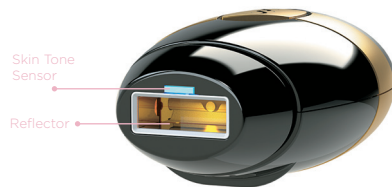10. TRAVELING WITH THE SMOOTHSKIN BARE+

The Power Supply unit of the SmoothSkin Bare+ is universal and will work with all commonly encountered mains supplies. Providing you have a suitable adaptor for the mains power cable plug, the SmoothSkin Bare+ will operate as normal.

11. TROUBLESHOOTING

PROBLEM	SOLUTION
The Activation Button Lights on the Handset do not light up when placed against your skin.	Make sure the power supply is plugged in. Ensure the Skin Tone Sensor is fully pressed against the skin. Your skin could be too dark for treatment. Check against the skin tone chart in the front of this user guide.
The Activation Button on the Handset lights up white when placed on the skin, but the Handset fails to flash when the Activation Button is pressed.	There could be a fault with the device. Unplug the device and restart. If the problem persists, visit the website for support details.
The Activation Button Lights on the Handset are illuminated red continuously, and the device will not function.	This indicates that there is a fault with the device. Stop using the device and visit the website for support details.
The SmoothSkin Bare+ Handset, Treatment Window, Power Supply or cable are damaged, broken, cracked or appear defective.	DO NOT USE. If you are in any doubt about the safety of SmoothSkin Bare+ or suspect it is damaged in any way, it must not be used, and you should visit the SmoothSkin website for more information.

12. WARRANTY

We grant a 2 year warranty on the product starting on the date of purchase. Within the warranty period, we will check faults in materials or workmanship free of charge and replace the complete device, at our discretion.

This warranty extends to every country where this device is supplied by CyDen or its appointed distributor.

This warranty does not cover damage due to improper use, normal wear or use as well as defects that have a negligible effect on the value or operation of the appliance. The warranty becomes void if repairs are undertaken by unauthorised persons and if original SmoothSkin Bare+ parts are not used. To obtain service within the warranty period, please visit our website: www.smoothskin.com.

To obtain service within the warranty period, please visit our website: www.smoothskin.com.

13. LABELS AND SYMBOLS

SmoothSkin Bare+ is marked with the following symbols:

	Type BF Applied Part		FCC EMC compliance for North America
	Class II Equipment		Australia Regulatory Compliance Mark
	Follow Instructions for use		EAC Conformity Marking for Custom Union
	WEEE waste disposal mark of the European Union		Trefoil compliance mark for Ukraine
	European Conformity Marking		China RoHS environmental protection use period
	Keep Dry		Warning Optical Radiation from device may cause eye injury
	NRTL Certification marking		

FCC Statement:

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

14. TECHNICAL SPECIFICATIONS

The SmoothSkin Bare+ is a filtered broadband Intense Pulsed Light system with the following technical specifications:

- **Repetition Rate:** Manually pulsed every 0.6 seconds. Continuous operation
- **Max Optical Output:** 4J/cm²
- **Pulse Length:** 1msec FWHM
- **Wavelength:** 510-1100nm
- **Power Line Input:** 100-240V~, 50/60Hz, 1.7-0.91A
- **Operating temperature:** Between 5°C and 40°C
- **Operating humidity:** Up to 93%R.H. non condensing
- **Operating Pressures:** 700hPa to 1060hPa
- **Treatment area (spot size):** 3 cm² (30mm x 10mm)

The area directly around the treatment window is classified as the Applied Part, and can reach a maximum temperature of 52degC.

Ensure SmoothSkin Bare+ is used as directed in Section 8.

Transport Conditions:

- **Temperature:** -25°C to +70°C
- **Humidity:** Up to 93%R.H. non condensing
- **Pressures:** 500hPa to 1060hPa

If the device has been stored at temperatures outside of the stated operating range, then please allow at least 30 minutes before use to allow the device to acclimatise.

SmoothSkin Bare+ has unlimited flashes and is intended for a single user. Minimum expected service life is 10 years when following the recommended regimen.

EMC

SmoothSkin Bare+ complies with the latest EMC Standards and Amendments applicable to Household and Medical Devices. These are: EN 60601-1-2:2014 Ed 4.0, EN 55014-1:2017 And EN 55014-2:2015, EN 61000-3-2:2014 & EN 61000-3-3:2013, FCC CRF 47 Part 15B 15.107 & 15.109"

Safety

SmoothSkin Bare+ complies with the latest Safety Standards and Amendments applicable to Household and Medical Devices.

These are: IEC 60335-1:2010, IEC 60335-2-23 :2003, IEC 60601-1:2015 , IEC 60601-1-1:2005, IEC 60601-2-57:2011, IEC 60601-1-6:2010, IEC 62366-1:2015, IEC 62366-1:2015, IEC 62304:2006, CAN/CSA -C22.2 No 60601-1:2008 and 2014, ANSI/AAMI ES 60601-1:2005 + AMS"

CE marking certifies that this appliance conforms to the following EEC directives:

- Low Voltage Directive 2014/35/EU
- Electromagnetic Compatibility Directive 2014/30/EU
- Energy related Product Directive 2009 / 125 / EC

DISPOSAL (end of life)

In order to minimize hazards to health and the environment and ensure that materials can be recycled, this product should be disposed of at a separate collection facility for waste electrical and electronic equipment.

The WEEE symbol marked on the product is there to remind you.

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Si desea más información, póngase en contacto con CyDen Ltd
CyDen Ltd,